



22 mm Flextube™ breathing system with integral monitoring line



Anaesthesia ▪ Anaesthetic Breathing Systems

22 mm Flextube™ breathing system with integral monitoring line

This range of 22 mm breathing systems are supplied with a gas monitoring line conveniently and safely positioned inside the length of the breathing system, providing a safer option.

By anchoring the monitoring line securely inside the breathing system optimal positioning is achieved for accurate and reliable gas monitoring and the need for additional luer ports in either the Y-Piece or filter/HMEF device is eliminated. This therefore avoids any possibility of accidental disconnection in use or misconnections.

The monitoring line is kept clear of the working environment, providing a tidy, less cluttered anaesthetic space and avoids the risk of occlusion or disconnection.

Combining the two products provides the convenience of a breathing system set up and ready to use and benefits the clinician by:

- providing consistent and accurate monitoring
- reducing the number of patient connections and the associated risks
- declutters the anaesthetic space around the patient

The Flextube breathing system with integral monitoring line is available in two lengths and two different configurations to suit clinical requirements.


Convenience

One product, complete and ready to use

Accuracy

Ensures consistent and accurate CO₂ monitoring

Optimal positioning of the monitoring line at the Y-piece

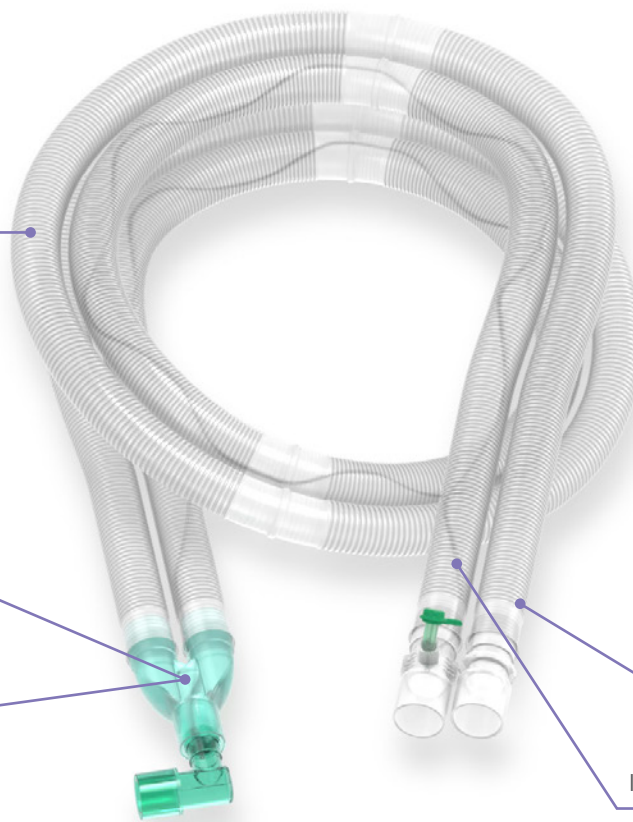

Safety

Helps to reduce the number of patient connections

Minimises risk of disconnection and occlusion



Safe and Secure

Flextube breathing systems have a unique clip feature to create an airtight fit when used with Intersurgical connectors, ensuring no leaks or disconnections

Reduces possibility of monitoring line becoming caught or entangled

Code	Description	Box Qty.
2010257	22 mm Flextube breathing system with integral monitoring line, ≥1.6 m	12
2010258	22 mm Flextube breathing system with integral monitoring line, ≥2.4 m	8
2010259	22 mm Flextube breathing system with 2 L bag, integral monitoring line and limb, ≥1.6 m	10
2010260	22 mm Flextube breathing system with 2 L bag, integral monitoring line and limb, ≥2.4 m	8

Make an enquiry [▶ Watch the video](#)

IS5.34 ▪ Issue 2 08.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.