



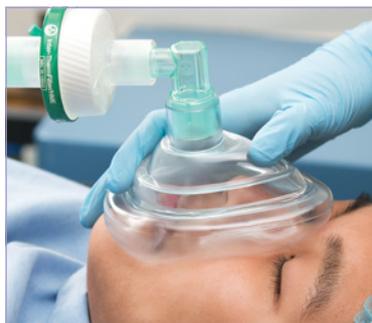

Anaesthetic face masks without a hook ring

A clear choice to reduce your impact on the environment

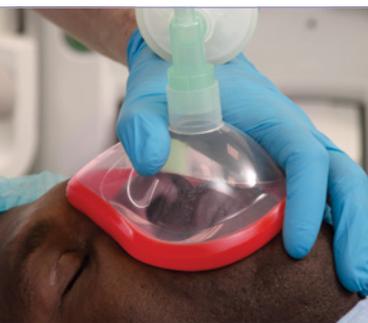

Anaesthesia ▪ Anaesthetic Face Masks



ClearLite™ anaesthetic face masks

A range of anaesthetic face masks which includes seven sizes to meet the clinical needs of all patient types and sizes. The range is available in sizes 0 neonatal to size 6 large adult. Due to the enhanced anatomical design, the ClearLite mask offers superior seal performance with minimum pressure, reducing the risk of unnecessary discomfort or potential trauma. The superior seal will also assist the clinician in

providing optimum manual ventilation of the patient in both the anaesthetic and emergency situation.

ClearLite is non-PVC and phthalate free which significantly minimises the impact upon the environment.

Lightweight mask

Provides enhanced grip and comfort in use for the clinician

Anatomical shape

Provides a precise fit for different face types and sizes. Minimum pressure required to achieve an effective seal reducing the risk of inadvertent trauma and nerve damage

Crystal clear body

Provides excellent visibility for confirmation that the patient is well perfused

Flexible non-inflatable cuff

Provides a comfortable and effective airtight seal



Safe and secure connection

This is achieved between the mask and the breathing system in use. The ISO 22mm tapered connection complies with ISO 5356

Reduced mask volume

Reduces potential rebreathing and risks associated with elevated levels of end tidal CO₂

Colour-coding

The comfortable, flexible seal is colour coded to allow quick and easy identification of the correct sized mask

Low environmental impact

Non-PVC and phthalate free

Video available at www.intersurgical.com



| Code | Description | Box Qty. |
|---------|---|----------|
| 7296001 | ClearLite, anaesthetic face mask, size 6, extra large adult, red seal, 22F | 20 |
| 7295001 | ClearLite, anaesthetic face mask, size 5, large adult, orange seal, 22F | 35 |
| 7294001 | ClearLite, anaesthetic face mask, size 4, adult, green seal, 22F | 35 |
| 7293001 | ClearLite, anaesthetic face mask, size 3, small adult, yellow seal, 22F | 35 |
| 7292001 | ClearLite, anaesthetic face mask, size 2, paediatric, white seal, 22F | 25 |
| 7291000 | ClearLite, anaesthetic face mask, size 1, infant, grey seal, no hook ring, 15M | 30 |
| 7290000 | ClearLite, anaesthetic face mask, size 0, neonate, blue seal, no hook ring, 15M | 40 |

[Make an enquiry](#)



Economy anaesthetic face masks

A range of seven colour-coded economy PVC anaesthetic face masks with a soft cushion and clear housing. Seven sizes from neonatal to extra large adult.



22F

Clear flexible housing and cushion

For good patient visibility

Soft cushion

For patient comfort and air tight fit

Colour coded packaging

for quick and easy identification

Seven sizes available

From neonatal to extra large adult



| Code | Description | Box Qty. |
|-----------|--|----------|
| ● 1517111 | Economy, anaesthetic face mask, size 6, extra large adult, 22F | 25 |
| ● 1516111 | Economy, anaesthetic face mask, size 5, large adult, 22F | 30 |
| ● 1515111 | Economy, anaesthetic face mask, size 4, medium adult, 22F | 40 |
| ● 1514111 | Economy, anaesthetic face mask, size 3, small adult, 22F | 50 |
| ○ 1513111 | Economy, anaesthetic face mask, size 2, paediatric, 22F | 35 |
| ● 1512111 | Economy, anaesthetic face mask, size 1, infant, 15M | 40 |
| ● 1511111 | Economy, anaesthetic face mask, size 0, neonate, 15M | 50 |

Make an enquiry

Scented anaesthetic face masks

In a range of sizes and scents

A range of scented single use PVC anaesthetic face masks with adjustable cushion, available in seven sizes and three scents: cherry, strawberry or vanilla, presented in a cool mint colour to help the patient relax.



Cushion valve at top (nose) of face mask

For user comfort

Soft adjustable face cushion

For patient comfort and a good anatomical fit



Colour coded packaging
for quick and easy identification

Clear flexible housing
For good patient facial visibility

Scented
To help patient relax

| Code | Description | Box Qty. |
|-----------|---|----------|
| ● 1126111 | Scented, vanilla, anaesthetic face mask, size 6, extra large adult, 22F | 24 |
| ● 1125111 | Scented, vanilla, anaesthetic face mask, size 5, large adult, 22F | 30 |
| ● 1124111 | Scented, vanilla, anaesthetic face mask, size 4, medium adult, 22F | 40 |
| ● 1123111 | Scented, vanilla, anaesthetic face mask, size 3, small adult, 22F | 50 |
| ○ 1122111 | Scented, cherry, anaesthetic face mask, size 2, paediatric, 22F | 35 |
| ○ 1129111 | Scented, strawberry, anaesthetic face mask, size 2, paediatric, 22F | 35 |
| ● 1121111 | Scented, cherry, anaesthetic face mask, size 1, infant, 15M | 32 |
| ● 1128111 | Scented, strawberry, anaesthetic face mask, size 1, infant, 15M | 32 |
| ● 1120111 | Scented, cherry, anaesthetic face mask, size 0, neonate, 15M | 40 |
| ● 1127111 | Scented, strawberry, anaesthetic face mask, size 0, neonate, 15M | 40 |

Make an enquiry



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.

IS5.35 ▪ Issue 1 02.20