



Intersurgical EcoLite™ paediatric mask range



Oxygen and Aerosol Therapy • Paediatric Masks



Intersurgical EcoLite™ paediatric mask range

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the mask's shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patient's face.

Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask body and seal, reducing the environmental impact.

Intersurgical EcoLite

[Find out more](#)

Eco range

[Find out more](#)



Incurved nose seal


Improving the fit and eliminating the need for a separate nose clip. Designed to minimise the risk of oxygen entering the patient's eyes

Soft face seal

For improved patient comfort

A choice of "under chin" or "on-chin" positions

Designed to fit on a wide range of patient faces

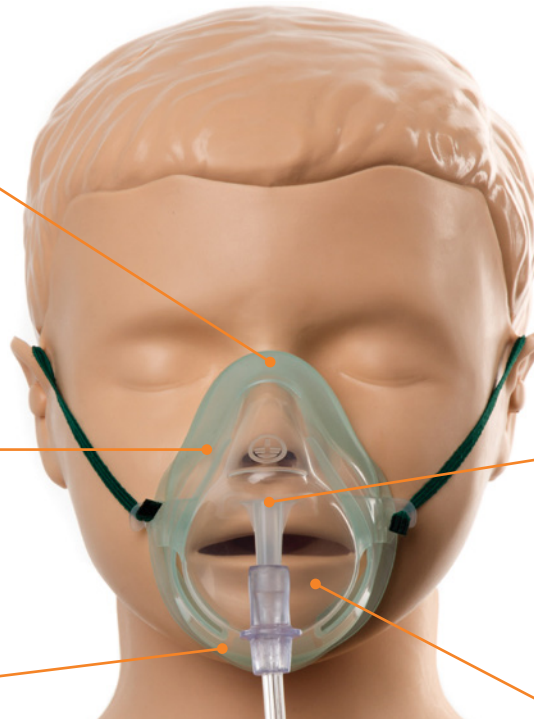

Lighter

Mask made with lightweight materials

Greener

The elimination of PVC from the masks has lowered their environmental impact

The shape and position of the exhalation ports have been designed to minimise the amount of rebreathed carbon dioxide

The material forming the body of the mask is clear, permitting visualisation of the patient's nose and lips. It is also rigid enough to maintain the mask's shape during transportation and use



Available with oxygen tubing featuring a universal connector, which can be attached to a 6mm flow meter spigot or fitted directly onto a 9/16-18 UNF Diameter Index Safety System (D.I.S.S.) flow meter



Paediatric medium concentration oxygen mask

The Intersurgical EcoLite paediatric medium concentration oxygen mask is extremely lightweight, and comfortable.

[Make an enquiry](#)
[Watch the video](#)


Code	Description	Tube length	Box Qty.	
1190015	Intersurgical EcoLite, paediatric, medium concentration oxygen mask		42	
1196015	Intersurgical EcoLite, paediatric, medium concentration oxygen mask with tube	2.1 m	40	
1196010	Intersurgical EcoLite, paediatric, medium concentration oxygen mask with non-PVC tube	2.1 m	40	
1196019	Intersurgical EcoLite, paediatric, medium concentration oxygen mask and tube with universal connector	2.1 m	40	

Paediatric adjustable venturi mask kits

The Intersurgical EcoLite™ paediatric adjustable venturi mask kits, available with the Silente™ and MultiOx™ adjustable venturi valve.

[Make an enquiry](#)


Code	Description	Tube length	Box Qty.	
1138015	Silente, Intersurgical EcoLite, paediatric, adjustable venturi mask kit with tube	1.8 m	40	
1118005	MultiOx, Intersurgical EcoLite, paediatric, adjustable venturi mask kit with tube	1.8 m	40	







Paediatric aerosol face masks and nebulisers

The Intersurgical EcoLite paediatric aerosol mask is available individually or with Cirrus™2 or HOT Top™2 nebuliser.

Make an enquiry




Code	Description	Tube length	Box Qty.
1198015	Intersurgical EcoLite, paediatric, aerosol mask		42 
1454015	Cirrus2 nebuliser, paediatric, Intersurgical EcoLite mask kit and tube	2.1 m	36 
1454010	Cirrus2 nebuliser, paediatric, Intersurgical EcoLite mask kit with non-PVC tube	2.1 m	36 
1468015	HOT Top2 nebuliser, paediatric, Intersurgical EcoLite mask kit with tube	1.8 m	36 

Venturi valve Intersurgical EcoLite™ mask kit

This product incorporates the Intersurgical EcoLite mask, offering patient comfort with a reduced environmental impact.

Make an enquiry



Code	Description	Tube length	Box Qty.
1107012	Intersurgical EcoLite, paediatric, venturi mask kit	1.8 m	30 

 Lower environmental impact product

IS10.27 • Issue 13 09.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.